

**Maximum Performance • Minimum Gas**

EOGas 4 reaches a Sterility Assurance Level (SAL) of  $10^{-6}$  in 3 hours using only 17.5g of ethylene oxide (EtO).




► **EOGas 4 Compared to Traditional EtO chambers:**

EOGas uses *less gas*, only 17.5g per cycle. Traditional chambers can use over 5 lbs of EtO per cycle, fed by a large external tank.

Traditional chamber sterilizers mix and dilute EtO with *inert gases*. The EOGas system uses 100% ethylene oxide.

The EOGas system performs sterilization and aeration in the *same chamber*. No need to transfer products to a separate aeration area while they are off-gassing.

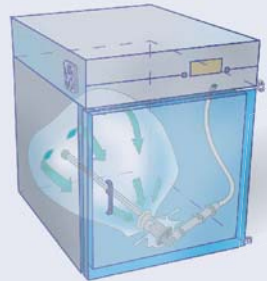

With EOGas sterilization and aeration occurs in the *same chamber*. There is no need to transfer products to a separate aeration area while gas is still being released.

► **EOGas 4 compared to other Low Temperature Sterilizers:**

The widest range of material and product compatibility (you can sterilize anything but foods, liquids, or drugs).


Ability to sterilize long, narrow, and multi-channel lumens.

All sides and surfaces are sterilized; No “shadowing”.

Compatible with all common packaging materials.


Will not damage delicate instruments like harsh chemical sterilants.

**EtO concentration in an EOGas 4 showing a SAL of  $10^{-6}$**

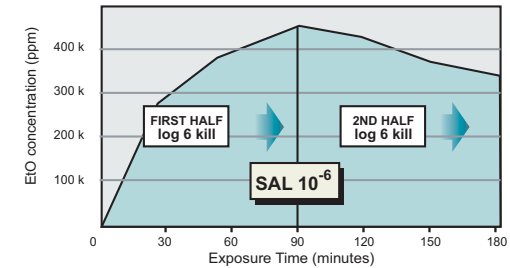


For more, see: [www.anpro.com/eogas4](http://www.anpro.com/eogas4)

**EOGas 4 Series vs. Sterrad: Maximum recommended length for sterilizing a 1mm Steel Lumen**

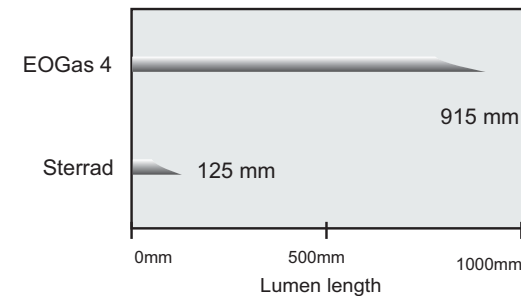

For our lumen study and comparison, see: [www.anpro.com/eogas4](http://www.anpro.com/eogas4)

► Specifications and Accessories

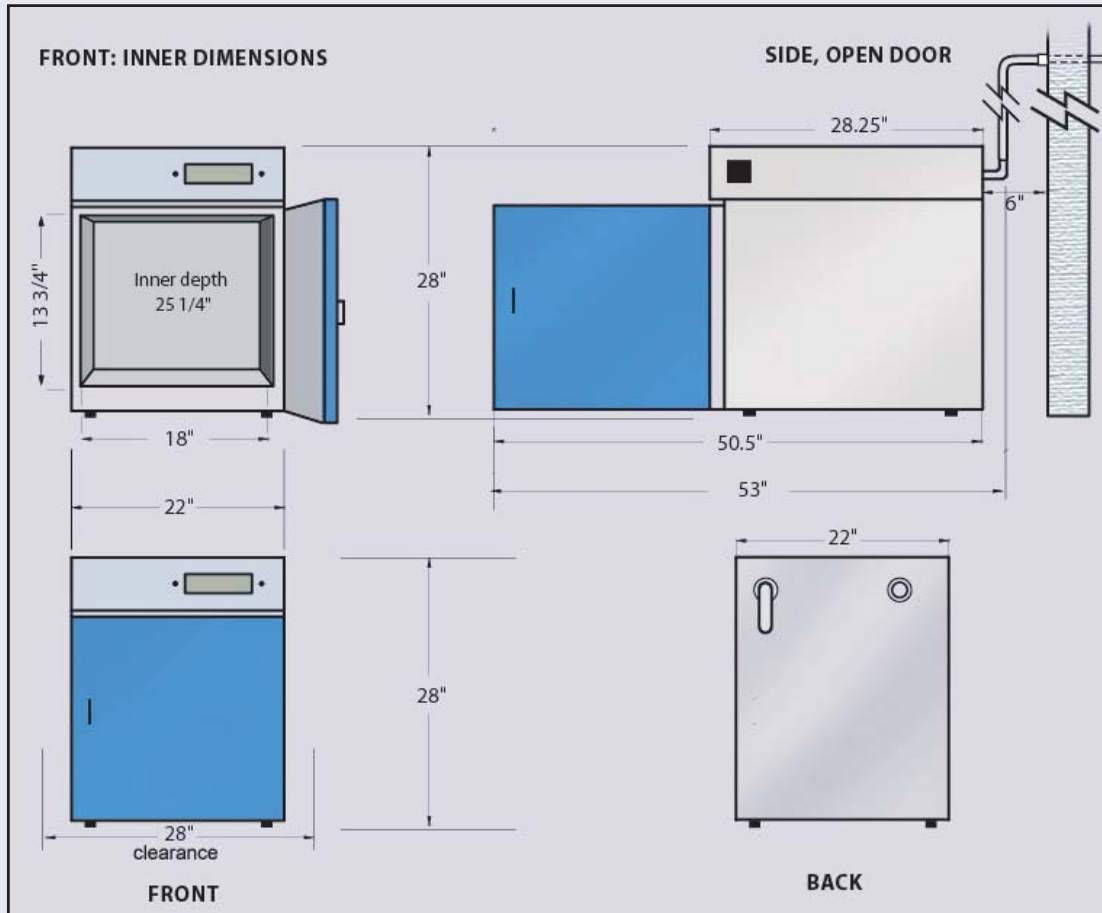

Model: AN4000.01, AN4000.11  
 Exterior dimensions: See diagram above  
 Inner dimensions: See diagram above  
 Total Weight: 168 lbs.  
 Power Supply: 115+/-10%, 150 Watts, 230+/-10%, 150 Watts  
 Manufactured by: Andersen Sterilizers, Inc.  
 Distributed by: Andersen Products, Inc.

**Starter Kit:**  
 Each Kit includes:  
 25 100% EtO cartridges  
 25 Humidichips  
 25 EOGas Dosimeters  
 25 EOGas bags

► Technical Support



**Operator Certification and Training**

Free training is provided for each EOGas operator. Training focuses on instrument preparation, packaging, operation and safety.



*NOTE: EOGas 4 has not yet been reviewed by the FDA for use in human medical facilities. Please contact us regarding any questions about products that may require extended aeration.*

Call us:  
**800-523-1276**  
 Or visit:  
**[www.anpro.com](http://www.anpro.com)**

**Andersen Products, Inc.**  
 Health Science Park  
 3202 Caroline Drive  
 Haw River, NC 27258-9564 USA  
 800-523-1276 • fax 336-376-8153

